



Tinzyme Co., Limited

Email: sales@tinzyme.com

Website: www.tinzyme.com

Tel: +86-755-86134126

WhatsApp/Facebook/Twitter: +86-189-22896756

FlashCut™ TaqI

Product Number: RE0574

Shipping and Storage

-20°C.

Components

Component	Specifications
FlashCut™ TaqI	200μL
10× FlashCut™ Buffer	2×1mL
10× FlashCut™ Color Buffer	2×1mL

Description

FlashCut™ Rapid endonucleases are a series of genetically engineered restriction endonucleases that are suitable for rapid enzymatic cleavage of plasmid DNA, PCR products, or genomic DNA. All FlashCut™ Rapid endonucleases have excellent activity in both FlashCut™ and FlashCut™ Color Buffer, and can complete enzyme cleavage within 5-15 minutes. In addition, Baishimei dephosphorylation and ligation reagents are available on FlashCut™ Buffer has 100% activity and supports one tube reaction, enhancing the experience of "enzyme digestion modification connection".

FlashCut™ Color Buffer includes red and yellow tracer dyes, which can be directly used for gel electrophoresis. The migration rate of red dye of FlashCut™ Color Buffer and 2500 bp double stranded DNA fragment in 1% agarose gel is similar; The migration rate of yellow dye and 10 bp double stranded DNA fragment in 1% agarose gel is similar.

Suggested reaction conditions

1. 1 × FlashCut™ Buffer solution.
2. Incubate at 37°C.
3. Prepare the reaction system according to the "DNA rapid enzyme digestion process".

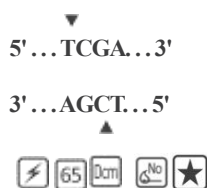
Inactivation conditions






Cannot be thermally deactivated, please use phenol chloroform extraction or column purification.

Quality control

1. Functional activity detection: At 65°C, in a 20μL universal FlashCut™ reaction system, 1μL FlashCut™ TaqI can completely digest 1μg λ DNA(Dam -) within 15 minutes.
2. Long term incubation detection: At 65°C, 1μL FlashCut™ TaqI was incubated with 1μg λ DNA (Dam -) in a 20μL universal FlashCut™ reaction system for 3 hours, and no non-specific degradation of the substrate caused by other nuclease contamination or star activity was detected. Longer enzyme digestion may result in star activity.
3. Enzyme digestion ligation re digestion detection: At 65°C, use FlashCut™ TaqI with 10 times the enzyme amount digests DNA substrates, recovers enzyme cleavage products, and at 22°C, T4 DNA Ligase (Fast) can reconnect over 95% of the enzyme cleavage products. After recycling the connecting product again, more than 95% of the connecting product can be re cut using the same endonuclease.

Icon annotation



1.  Rapid endonuclease can complete the reaction within 5-15 minutes.
2.  The optimal reaction temperature is 65°C.
3.  For DNA methylated by Dcm, splicing may be hindered
4.  Cannot be thermally deactivated.
5.  3 hours of incubation did not show star activity, and longer enzyme digestion may result in star activity.

Protocol

1. DNA rapid enzymatic digestion process

1.1. Prepare the reaction system on ice according to the recommended sample addition sequence as follows:

	plasmid DNA	PCR products	Genomic DNA
ddH ₂ O	15μL	16μL	30μL
10×FlashCut™ Buffer or 10×FlashCut™	2μL	3μL	5μL
Color Buffer			
Substrate DNA	2μL(up to 1μg)	10μL(up to~0.2μg)	10μL(5μg)
FlashCut™ TaqI	1μL	1μL	5μL
Total	20μL	30μL	50μL

This system is suitable for enzyme digestion of purified PCR products. The unpurified PCR product has a certain ionic strength, 10×FlashCut™. The amount of buffer added can be appropriately reduced to 2μL. However, due to the simultaneous exonuclease activity of DNA polymerase, it can affect the cleavage products. Therefore, the following steps require cloning and other operations. It is recommended to purify the PCR products before enzyme cleavage.

1.2. Gently suck or tap the tube wall to mix well (without vortexing), then centrifuge instantly to collect hanging wall droplets.

1.3. Incubate at 37°C for 15 minutes (plasmid), or 15-30 minutes (PCR product), or 30-60 minutes (genomic DNA).

Note: Enzyme activity may slightly decrease during the reaction at 37°C, and the incubation time should be appropriately extended.

1.4. Phenol chloroform extraction or column purification (optional).

1.5. If using FlashCut™ Color Buffer is used for enzyme digestion reaction, and the resulting product can be directly subjected to sample electrophoresis.

2. Double enzyme digestion or multi enzyme digestion

2.1. The dosage of each rapid endonuclease is 1μL, and the reaction system should be appropriately expanded as needed;

2.2. The total volume of all rapid endonucleases must not exceed 1/10 of the total reaction system;

2.3. If the optimal reaction temperatures for the several rapid endonucleases used are different, the enzyme with the lower optimal temperature should be used first for enzyme digestion, and then the enzyme with the higher optimal temperature should be added for enzyme digestion reaction at its optimal reaction temperature.

3. Expansion reaction system suitable for plasmids

DNA	1μg	2μg	3μg	4μg	5μg
FlashCut™ TaqI	1μL	2μL	3μL	4μL	5μL
10×FlashCut™ Buffer or 10×FlashCut™	1μL	2μL	3μL	4μL	5μL
Color Buffer					
Total	20μL	20μL	30μL	40μL	50μL



Tinzyme Co., Limited

Email: sales@tinzyme.com

Website: www.tinzyme.com

Tel: +86-755-86134126

WhatsApp/Facebook/Twitter: +86-189-22896756

Note: If the total reaction system is greater than 20 μ L, the incubation time should be appropriately increased, and water bath, metal bath or sand bath should be used as much as possible.

The number of enzyme cleavage sites in different DNA

λ DNA	Φ X174	pBR322	pUC57	pUC18/19	SV40	M13mp18/19	Adeno2
121	10	7	4	4	1	12	50

The influence of methylation modification

Dam	Dcm	CpG	EcoKI	EcoBI
Shear obstruction	no effect	no effect	no effect	no effect

Activity in different reaction buffers

	Cut One™ Buffer	Thermo Scientific Fast Digest Buffer	NEB Cut Smart® Buffer	Takara Quick Cut™ Buffer
Reactivity	100%	100%	100%	100%

Note: The activity data comes from the detection under the restriction enzyme standard reaction system.

The activity of DNA modifying enzymes in FlashCut™ Buffer and FlashCut™ Color Buffer

Alkaline Phosphatase (Fast)	100%
T4 DNA Ligase (Fast)	100%

Note: The activity data comes from the detection under standard reaction system, and T4 DNA Ligase (Fast) requires ATP as a cofactor.